

FluoroDose[®]

5% Sodium Fluoride Varnish with Xylitol



✓ **Easier to hold easier to use**

✓ **Tastes Better**

“Every child liked it!” and “Patients like the variety in flavors.” Flavors are tasty and not medicinal. Xylitol sweetener provides added anti-bacterial properties with no bitter aftertaste

✓ **Feels Better**

“Application is smoother than most varnishes.” Patients appreciate the almost invisible texture. FluoroDose dries in seconds and adheres to the teeth for up to six hours for maximum uptake

✓ **Dries clear, smooth**


“Applies smoothly without clumping and dries clear.” Patients leave with natural-looking teeth; no yellowing or discoloration. It’s less sticky and won’t clog evacuation systems

✓ **Perfect Size**

“Just the right amount for a single-patient treatment without waste.” Single-patient, 0.3 mL dose eliminates cross-contamination and provides just the right amount for kids or adults.

✓ **Best Value**

“93% Clinical Rating. 92% Recommend It. 81% Switched to It.” FluoroDose is the best value fluoride varnish available today.

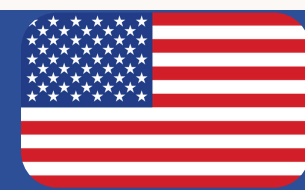

FluoroDose[®]

5% Sodium Fluoride Varnish with Xylitol

✓ Indications For Use

- Treats tooth sensitivity to heat, cold and sweets when applied to exposed roots.
- Helps prevent root decay in older patients who experience difficulty in cleaning their teeth by occluding dentin tubules.
- The varnish application should be repeated at three-month intervals for high-risk patients and at six-month intervals for patients who are not at high risk.

✓ Allergen Risk Information

FluoroDose does not contain lactose, nuts, dairy, eggs, gluten, soy, peanuts, tree nuts or their derivatives, fish/shellfish, seed products or any dyes (including Red Dye 40).

FluoroDose does contain rosin and beeswax

✓ Flavors

- Bubble Gum
- Caramel
- Cheery
- Melon

✓ Warning

- Keep away from heat, sparks, open flame and hot surfaces.
- Keep out of the reach of children.

✓ Packing

Each unit contains one single dose (0.3 mL) of FluoroDose plus one Benda[®] Brush applicator.

Simply peel back the cover, dip in the brush, mix and apply.

✓ **To receive the maximum benefits of treatment**, the varnish needs to remain on your teeth for the next 4 to 6 hours. During this time, we recommend the following: Patient should be advised to remain on a soft food diet and only drink cold liquids for 2 hours after treatment; and not to brush or floss for a minimum of 4 to 6 hours. Avoid alcoholic beverages.

PREVENTION FOR LIFE[®] = CARIES FREE FOR LIFE!

